

MicroFlow III Ductless Workstation

The MicroFlow III Workstation ideal for: A class 1 ductless hood workstation Equipped with activated carbon filtration, for fumes, odors, non-hazardous chemical vapors and particulate filtration.

- Completely self contained with integral recessed work surface to contain spills.
- Clear hood surround with safety viewing sash allows visibility from 3 sides, and provides both safety and access to inside work area. Sash can be conformed for use with a microscope.
- Variable speed fan control. High speed for 100f/m. Air flow thru sash opening.
- Dimensions 24" wide X 20.75" high X 24" deep.



Description	Cat. No.
1. MicroFlow III Workstation 115V	74010
2. Replacement carbon filter	52006
3. Replacement pre-filter	74011
4. Mobile Table	

- Operates on 115V AC, or 230V International, conforms to UL, CSA, & CE requirements.
- Electrical cord port exits both sides.
- Easy change out of used filter. Recommend filter change every 60-90 days. Filter life varies with usage.


Easy Filter Change



Remove plastic wrap from carbon filter & snap into place using filter clips



Install particle prefilter and attach using velcro

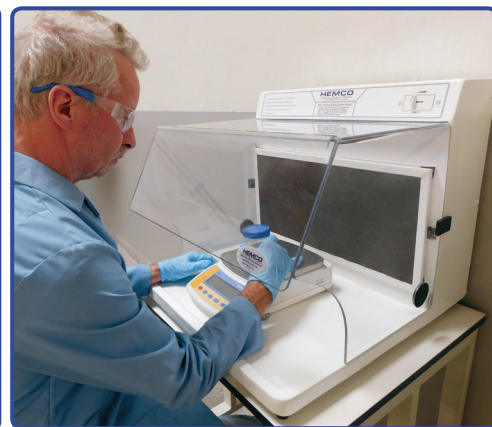

MicroFlow III is ready for use

HEMCO Corporation
711 S. Powell Road
Independence, MO 64056

HEMCO
Lab Planning Solutions

Phone (816) 796-2900
Fax (816) 796-3333
MFIII-04.18